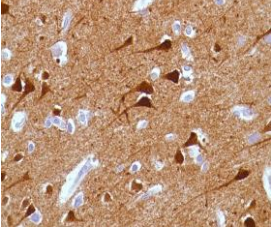




Rabbit Anti-Human Syntaxin 1A Monoclonal Antibody (Clone SP36)

CATALOG #:	M3360 0.1 ml rabbit monoclonal antibody purified by protein A/G in PBS/1% BSA buffer pH 7.6 with less than 0.1% sodium azide. M3362 0.5 ml rabbit monoclonal antibody purified by protein A/G in PBS/1% BSA buffer pH 7.6 with less than 0.1% sodium azide. M3364 1.0 ml rabbit monoclonal antibody purified by protein A/G in PBS/1% BSA buffer pH 7.6 with less than 0.1% sodium azide. M3361 7.0 ml pre-diluted rabbit monoclonal antibody purified by protein A/G in TBS/1% BSA buffer pH 7.6 with less than 0.1% sodium azide.	 <p><i>Human cerebral cortex stained with anti-syntaxin 1A antibody</i></p>
INTENDED USE:	For Research Use Only. Not for use in diagnostic procedures.	
CLONE:	SP36	
IMMUNOGEN:	Synthetic peptide derived from N-terminus of human syntaxin 1A.	
IG ISOTYPE:	Rabbit IgG	
EPITOPE:	Not determined	
MOLECULAR WEIGHT	33kDa	
SPECIES REACTIVITY:	Human and mouse (tested). (See www.springbio.com for information on species reactivity predicted by sequence homology.)	
DESCRIPTION:	Syntaxin 1A (STX1A) is potentially involved in docking of synaptic vesicles at presynaptic active zones and may play a critical role in neurotransmitter exocytosis. Haploinsufficiency of STX1A may be the cause of certain cardiovascular and musculoskeletal abnormalities observed in Williams-Beuren syndrome (WBS), a rare developmental disorder. It is a contiguous gene deletion syndrome involving genes from chromosome band 7q11.23.	
APPLICATIONS:	Immunohistochemistry (IHC)	
IHC PROCEDURE:	Specimen Preparation: Formalin-fixed, paraffin-embedded tissues are suitable for use with this primary antibody. Deparaffinization: Deparaffinize slides using xylene or xylene alternative and graded alcohols. Antibody Dilution: If using the concentrate format of this product, dilute the antibody 1:100. The dilutions are estimates; actual results may differ because of variability in methods and protocols. Antigen Retrieval: Boil tissue section in 10mM citrate buffer, pH 6.0 for 10 min followed by cooling at room temperature for 20 min. Primary Antibody Incubation: Incubate for 30 minutes at room temperature. Slide Washing: Slides must be washed in between steps. Rinse slides with PBS/0.05% Tween. Visualization: Detect the antibody as instructed by the instructions provided with the visualization system.	
POSITIVE CONTROL:	Cerebral cortex	
CELLULAR LOCALIZATION:	Membrane	
STORAGE & STABILITY	Store at 2-8°C. Do not freeze. The user must validate any other storage conditions. When properly stored, the reagent is stable to the date indicated on the label. Do not use the reagent beyond the expiration date. There are no definitive signs to indicate instability of this product; therefore, positive and negative controls should be tested simultaneously with unknown specimens. If unexpected results are observed which cannot be explained by variations in laboratory procedures and a problem with the reagent is suspected, contact Technical Support at spring.tech@ventana.roche.com .	
WARNINGS & PRECAUTIONS:	<ol style="list-style-type: none">1. Avoid contact of reagents with eyes and mucous membranes. If reagents come into contact with sensitive areas, wash with copious amounts of water.2. This product is harmful if swallowed.3. Consult local or state authorities with regard to recommended method of disposal.4. Avoid microbial contamination of reagents.	