



Rabbit Anti-Human MAP4K4 Monoclonal Antibody (Clone SP177)


CATALOG #:

M4770 0.1 ml rabbit monoclonal antibody purified by protein A/G in PBS/1% BSA buffer pH 7.6 with less than 0.1% sodium azide.


M4772 0.5 ml rabbit monoclonal antibody purified by protein A/G in PBS/1% BSA buffer pH 7.6 with less than 0.1% sodium azide.

M4774 1.0 ml rabbit monoclonal antibody purified by protein A/G in PBS/1% BSA buffer pH 7.6 with less than 0.1% sodium azide.

M4771 7.0 ml pre-diluted rabbit monoclonal antibody purified by protein A/G in TBS/1% BSA buffer pH 7.6 with less than 0.1% sodium azide.

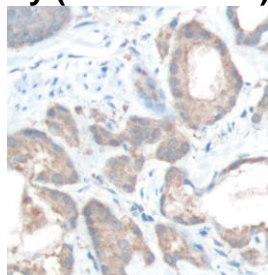

Human prostate adenocarcinoma stained with anti- MAP4K4 antibody

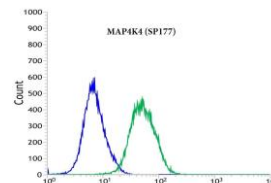

Flow cytometric analysis of rabbit anti-MAP4K4 (SP177) antibody in HEPG2 (green) compare to negative control of rabbit IgG (blue)

INTENDED USE:

For Research Use Only. Not for use in diagnostic procedures.

CLONE:

SP177

IMMUNOGEN:

Synthetic peptide derived from the C-terminus of human MAP4K4 protein.

IG ISOTYPE:

Rabbit IgG

EPITOPE:

Not determined

MOLECULAR WEIGHT:

142 kDa

SPECIES REACTIVITY:

Human (tested). (See www.springbio.com for information on species reactivity predicted by sequence homology.)

DESCRIPTION:

Mitogen-activated protein kinase kinase kinase 4 (MAP4K4) is a member of stress activated MAP kinase signaling pathway. It directly phosphorylates and activates the c-Jun N-terminal kinase (JNK) and p38 families of MAP kinases in response to various stimulating factors. It involves the mediation of TNF-alpha signaling pathway. MAP4K4 is overexpressed in many types of human cancer such as carcinoma of prostate, ovarian, liver, and pancreas.

APPLICATIONS:

Immunohistochemistry (IHC) and Flow Cytometry

IHC PROCEDURE:

Specimen Preparation: Formalin-fixed, paraffin-embedded tissues are suitable for use with this primary antibody.

Deparaffinization: Deparaffinize slides using xylene or xylene alternative and graded alcohols.

Antibody Dilution: If using the concentrate format of this product, dilute the antibody 1:100. The dilutions are estimates; actual results may differ because of variability in methods and protocols.

Antigen Retrieval: Boil tissue section in 10mM citrate buffer, pH 6.0 for 10 min followed by cooling at room temperature for 20 min.

Primary Antibody Incubation: Incubate for 10 minutes at room temperature.

Slide Washing: Slides must be washed in between steps. Rinse slides with PBS/0.05% Tween.

Visualization: Detect the antibody as instructed by the instructions provided with the visualization system.

IHC POSITIVE CONTROL:

Prostate Adenocarcinoma

FLOW CYTOMETRY:

Recommended starting protocol: Dilute the antibody 1:400. Incubate for 30 minutes at 4°C.

The dilution is an estimate; actual results may differ because of variability in methods and protocols. Optimal dilution and procedure should be determined by the end user.

FLOW CYTOMETRY

POSITIVE CONTROL:

HEPG2 Cell Line

CELLULAR LOCALIZATION:

Cytoplasm

STORAGE & STABILITY:

Store at 2-8°C. Do not freeze. The user must validate any other storage conditions. When properly stored, the reagent is stable to the date indicated on the label. Do not use the reagent beyond the expiration date.

There are no definitive signs to indicate instability of this product; therefore, positive and negative controls should be tested simultaneously with unknown specimens.

If unexpected results are observed which cannot be explained by variations in laboratory procedures and a problem with the reagent is suspected, contact Technical Support at spring.tech@ventana.roche.com.

**WARNINGS &
PRECAUTIONS:**

1. Avoid contact of reagents with eyes and mucous membranes. If reagents come into contact with sensitive areas, wash with copious amounts of water.
2. This product is harmful if swallowed.
3. Consult local or state authorities with regard to recommended method of disposal.
4. Avoid microbial contamination of reagents.